

The AccuPen[®]


Handheld Applanation Tonometer

When The Pressure's On...
The AccuPen Delivers.



Call to order the AccuPen today and receive 100 free AccuTip covers.

Product #
24-3000

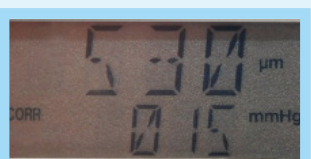

Busy facilities deserve an instrument that will come through in the clutch. The AccuPen's combination of accuracy and portability will make it an invaluable instrument to your entire staff.

Accurate

- Gravity Offset Technology provides precise IOP measurements with less calibration compared to other handheld tonometers.
- Micro strain gauge technology combined with proprietary algorithm provides reproducible results.

Easy-to-use

- Ergonomic design combined with sighting lines allow for easy visualization of the cornea, facilitating both centralization and perpendicularity.
- Contoured molding with large, visible LCD allows the AccuPen to be used ambidextrously and fit comfortably in any size hand.

Ready-to-use

- Eliminates the constant need for "up, down and around" calibration.
- Long lasting lithium battery lasts up to 15,000 readings.

Portability

- Weighs only 3.0 oz.
- Comes with custom, heavy duty traveling case and lanyard for safe transportation.

Versatility

- Takes fast measurements with patient in any position.
- Frees up valuable office space.

Smarter Software

- Stores and averages nine measurements.
- Built in IOP Correction Calculator adjusts IOP based on corneal thickness.

Specifications:	
Accuracy*	+/- 2mmHg
Weight	- 3.0 oz
Dimensions	- 7 1/4" x 1 1/4"
Transducer Dimensions	1.0 mm diameter
Measurement Range	5 - 60mmHg
Power Source	3.6 V Lithium Battery
*Compared to Goldmann Tonometer	